


CCP 3100

SYSTEM NETWORKING

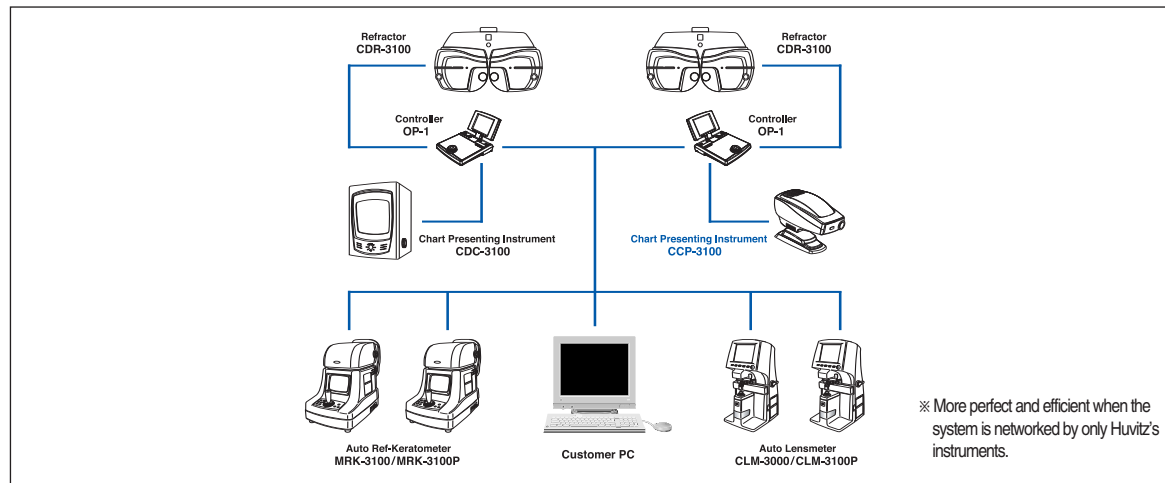

Basic Accessories

- Remote Control
- Screen (400X500mm)
- Table Stand
- Dust Cover
- Fuse
- Key Spanner
- Power Cable



SPECIFICATIONS

Chart	41 chart 34 mask, Red / Green & Polarization Filters
Projection Distance	2.5 ~ 8m
Chart Rotation Speed	Average 0.15 sec
Projection Magnification	30x at 5m
Power Saving	Automatic Switch off (10 min)
Program	2 programs with a maximum of 30 charts each
Tilt Angle	15°
Power Supply	110-120 / 200-240V ~, 50 / 60Hz, 0.6A
Lamp	6V / 20W (Halogen)
Dimension	367(W) X 219 (D) X 273 (H) mm
Weight	6.6 kg

Designs and details can be changed without prior notice for its improvements.

Huvitz

Huvitz Building 689-3, Geumjeong-dong
Gunpo-si, Gyeonggi-do, 435-862, Korea
Tel : ++82-31-442-8868
Fax : ++82-31-477-8617
http : //www.huvitz.com

Distributed by

CE

High quality Chart Projector which you want.

CCP 3100 CHART PROJECTOR



Huvitz
Pacing Progress toward People



The Huvitz CCP-3100 charts with high quality, allow precise acuity test for you.

The CCP-3100 delivers its promises by utilizing the total refraction system to perform any test you want accurately and efficiently.

The various kind of chart of the CCP-3100 can be changed more precisely within very short time.

The CCP-3100 satisfies all methods of optometry tests providing the best environment.

Various charts of the CCP-3100, You can choose any test you want.

High Quality Charts
Clear and Fast !

The CCP-3100 has high resolution power of line, 50lines/mm, so provides neat and clean graphics and a fast chart changing speed marking it easy for operators to use and comfortable for customers.


Wide range of features
Various and Great!

The CCP-3100 with 41 charts including Red/Green and Polarizing filters are the new solution for Binocular balance test, Stereo test, Aniseikonia test, Fusion test.


Great networking
Easy and Convenient !

Any refraction system including the Huvitz CDR-3100 can incorporate CCP-3100 to provide fast and efficient vision tests.


RC - A TYPE

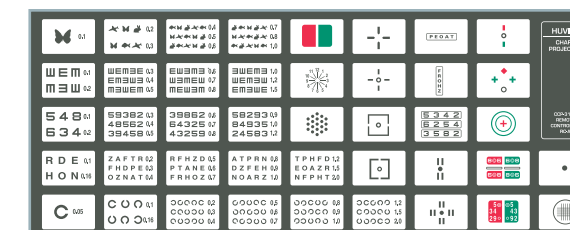

RC - C TYPE

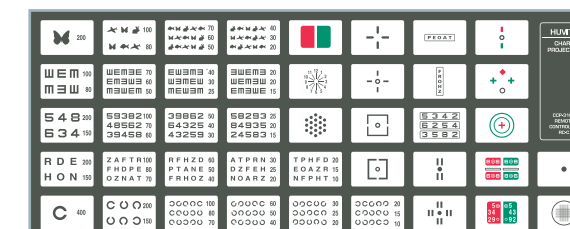

RC - D TYPE

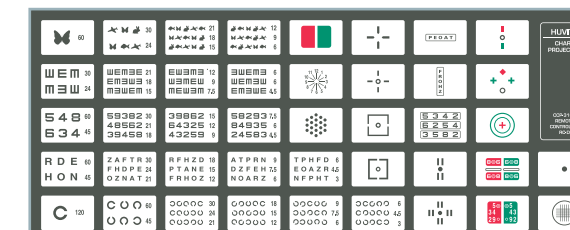

An elegant and compact Design

The CCP-3100 can be accommodated in small area. It can be stood up, hung on a wall, placed on a table stand.

Controllable by the CDR-3100 operation panel(OP-1)

Huvitz's high networking technology provides the remote control of CCP-3100 by the CDR-3100 operation panel(OP-1).